

# Solvent Safety Alarm System




## Solvent Safety Alarm System


- Monitor up to 3 containers (**Waste and/or Mobile Phase**) at the same time
- Sensors can be set to detect **FULL or EMPTY** fluid levels
- Front panel LEDs indicate which sensor position has been activated
- Beacon light illuminates to signify fluid level status:  
**Green**=Safe, **Red**=Alert!, **Orange**=Alert with audio alarm paused\*
- Adjustable audio alert




A sensor can be set to detect low fluid level in your mobile phase container.

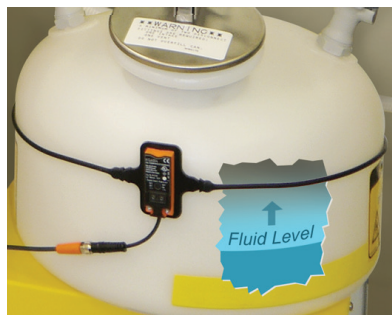

Another sensor can be set to detect high fluid level in your waste container.

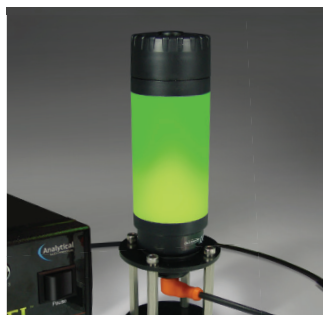

The Beacon turns from **green** to **red** and an audio alert sounds when a level has been reached (full or empty, depending on setup).

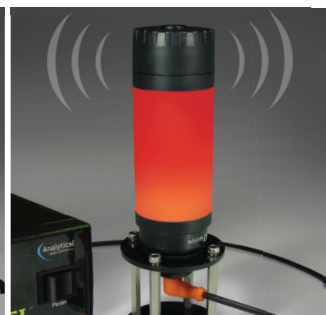

Cat. No.

Description

Sentinel

Solvent Safety Alarm System, Includes Control Module, Beacon and 3 Sensors



# Operation



A sensor monitors the fluid level in your container. The Beacon illuminates **green** under safe level conditions. Up to three containers can be monitored at once.




The Beacon illuminates **red** and an audio alert sounds when a level has been reached (full or empty, depending on setup).



Front panel LEDs indicate which sensor has been activated.



The **Pause** button feature suspends audio and visual (red alert) while the Beacon remains illuminated **orange**, reminding the user that attention is needed.

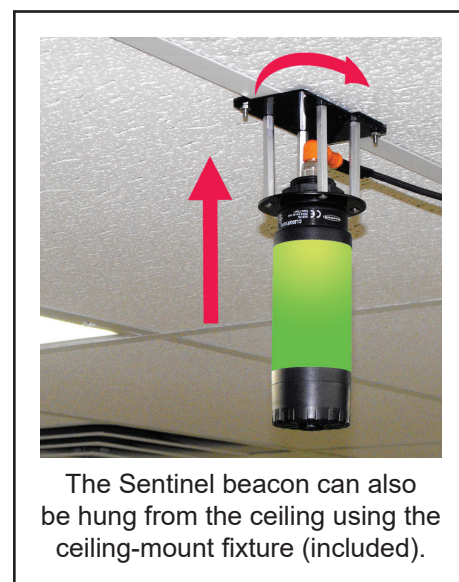

The Sentinel beacon can also be hung from the ceiling using the ceiling-mount fixture (included).