


X-Type Syringes

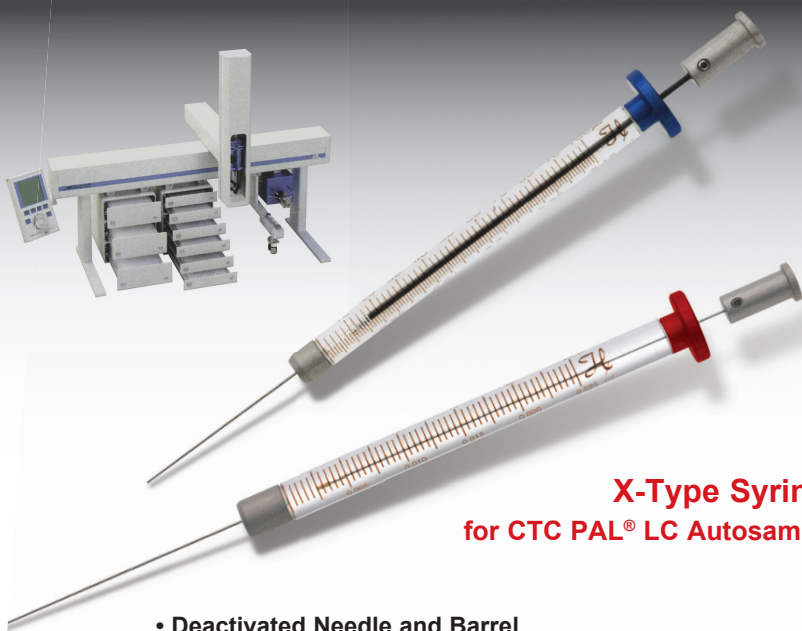

X-Type Syringes for CTC PAL® LC Autosamplers

- **Deactivated Needle and Barrel**

- Eliminates contact between samples and metal components
- Prevents sample adsorption on inner surface of the barrel
- Extremely smooth surface for low friction and increased plunger lifetime

- **Inert PTFE-Polymer Plunger Tip**

- High mechanical stability for increased syringe lifetime
- Highly inert plunger tip material compatible with organics and high and low pH solutions

- **Unique Needle Design**

- Needle attached directly to the barrel
- No glue or cement in fluid path
- Improved design for limited carryover

Fixed Needle (FN) Syringe (point style 3)

Cat. No.	Volume	Model	Gauge	CTC Part No.
42044754	25 µL	1702	22s	SYRX G25-22S-3
42044004	100 µL	1710	22	SYRX G100-22-3
42044524	100 µL	1710	22s	SYRX G100-22S-3

For pricing, see analytical-sales.com